

EMPRESA

— FEATURING FLOORBED TECHNOLOGY —

THE
EMINENCE
COLLECTION
by Accora



“Every design
begins with
an even
better story”

Lorinda Mamo, designer

Your story

At last, a bed which is a true all-rounder, combining your clinical, operational and design must-haves.

Your choice of nursing bed doesn't have to be a compromise. Even with the diverse demands of clinical care, staff and resident safety, robustness and stylish design – we've managed to bring these all together in our new **Empresa™** bed.

A bed for all.





Nursing beds transformed

First impressions count when it comes to attracting new residents. Looks play a huge part in making a home feel 'right' and just because specialist beds might be needed, it doesn't mean they can't be beautiful. With lots of options to choose from, our Empresa bed is a perfect partner for your design.

Working with our highly experienced team, we'll craft a centrepiece for you to be proud of – and help you create gorgeous spaces for people to live in.


Be inspired.

Exceptional safety

Accora has been helping make care safer since its inception. Our class beating FloorBed® technology brings proven effectiveness in bed falls and injury prevention¹, re-defining expectations of a floor level bed.

With falls as the leading cause of death and serious injury in older people², FloorBed technology can play an important role in falls management in longterm care settings. And where residents need extra protection, the addition of a high safety mat can help reduce risk even further, by transforming a potential fall into a simple roll.

Falls within the bedroom were reduced by 80% and falls from the bed were eliminated.¹

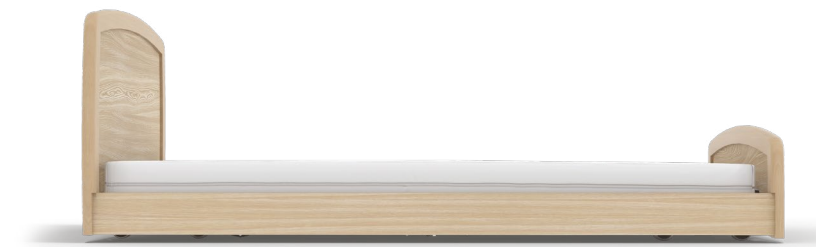

Lowers to 100 mm




High safety mat


Summary of product evaluation results

Description	Baseline	Post FloorBed	Trend
Residents who fell	9 (of 12)	5 (of 12)	▼ 44%
Falls in bedroom space (exc. bath)	15	3	▼ 80%
Notable injury and/or required diagnostic/medical intervention	2	0	▼ 100%
Falls from the bed	7	0	▼ 100%




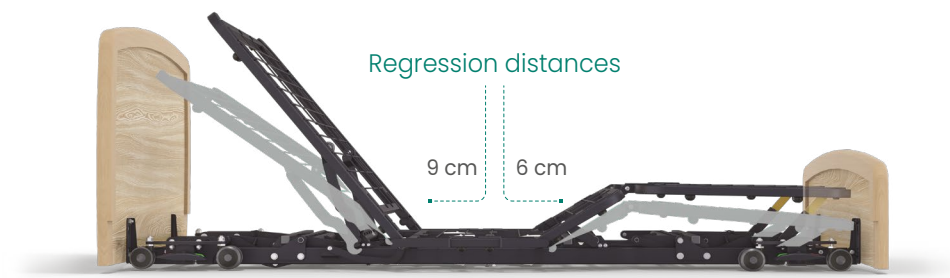


BodyMove regression distance

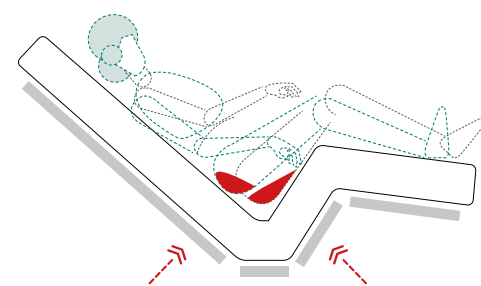


No regression increases friction and shear

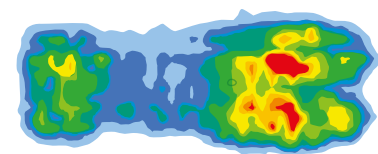

Pressure map showing no regression

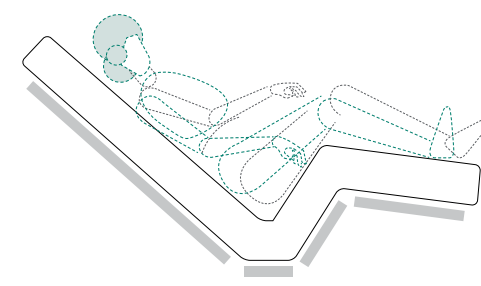

BodyMove regression reduces friction and shear

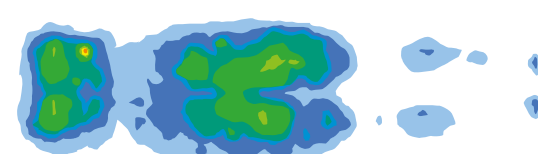

Pressure map showing BodyMove regression

Exceptional care

Caregiving is now complemented by the introduction of our dynamic **BodyMove™** system, helping form part of your first line of defence against tissue damage.

When moving from a supine to sitting position in bed, the body naturally lengthens. **Empresa's BodyMove** auto-regression feature accommodates this change, so not only does it help with pressure distribution, it helps improve comfort too. With all essential clinical features to aid safe nursing and residents' independence, the **Empresa** bed has people at the heart of its design.

Test results showed a 23% reduction in average pressure over the sacral area with the **Empresa BodyMove** technology, when compared to a standard nursing bed.³

“Good design can make the world of difference to how a resident, carer or relative will feel...”⁴

Jacqui Smith, interior designer



Ergonomic handling

A comfortable to hold handset with a simple one-touch operation, activates pre-programmed care and comfort positions, as well as other profiling options.

And with safety top of mind too, there is a lockout function to prevent unauthorised and inadvertent use, and an auto safety stop, inhibiting the bed from being lowered to below 200 mm. A double-button press is needed to lower the bed to the floor.

Helping ensure effortless manoeuvrability and even greater levels of safety, the **Empresa** bed also comes equipped with eight castors and our **UniSafe™** braking system. Moving the bed around the room couldn't be easier and it takes only a single press of a pedal to secure each set of brakes.

Putting you in control.



Care position



Comfort position



UniSafe™ braking system



Handset





Full nursing height 800 mm



1050 mm width adjustment



Built-in length extension

Ultimate flexibility

Everything you need in just one bed. Built-in adjustability gives the option to vary dimensions and a wide range of universal accessories allows you to utilise your assets fully, even when the profile of your resident changes. Combined with less risk of downtime, streamlined procurement and reduced training needs, the [Empresa](#) bed deserves a closer look.

And clinical functionality supports nursing care when needed, with full height adjustment for caregiver safety – all packaged in a choice of eye-catching and bespoke designs to help showcase your home.

Making life simpler.

“Details are
not details;
they make the
product”

Charles Eames, designer



Wood designs



Arc



Atelier



Bento



Chateau



Classic



Skandi



Vogue

Upholstered designs



Alexander



Alexander
with wood



Belgrave



Cleveland



Finsbury




Grosvenor

Wood finishes




Vicenza Oak H3157

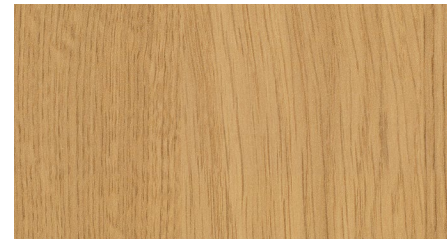


Natural Lancaster Oak H3368




Lissa Oak R20034

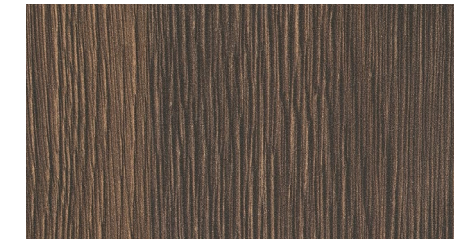


Tobacco Gladstone Oak H3325

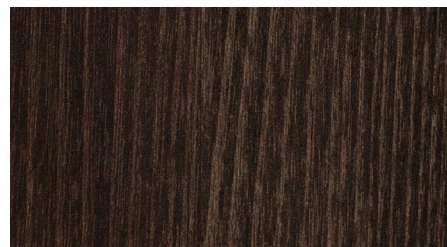


Black-Brown Thermo Oak H1199

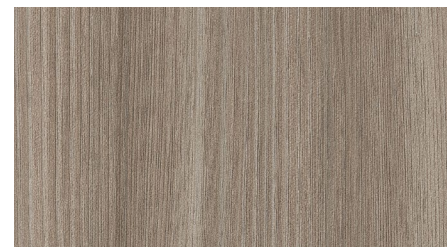

Grey-Beige Tossini Elm H1210

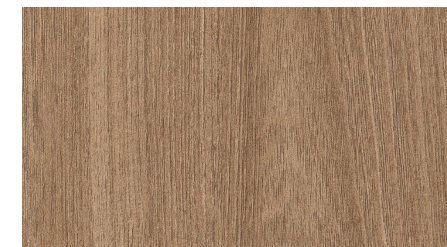

Brown Tossini Elm H1212

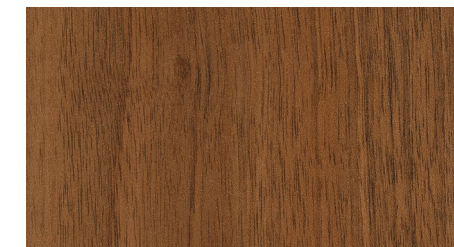


Natural Dijon Walnut H3734




Standard Walnut R30001




Lincoln Walnut H1714

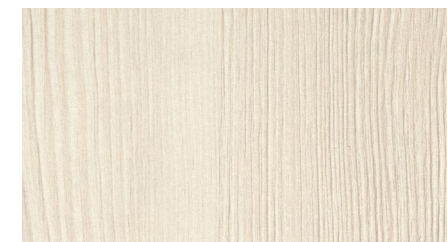


Polar Aland Pine H3433

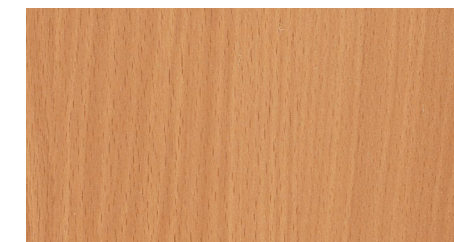


Bavarian Beech H1511

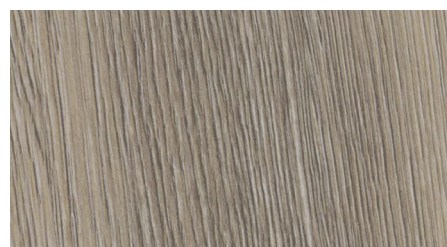


Grey Sacramento Oak H1143

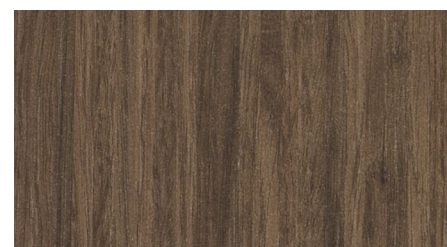

Brown Tonsberg Oak H309

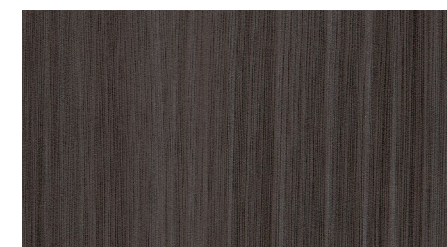


Anthracite Finline Metallic H3190

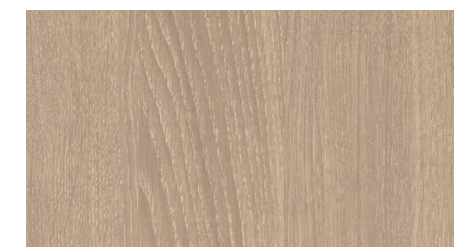


Sand Orleans Oak H1377

Upholstery finishes

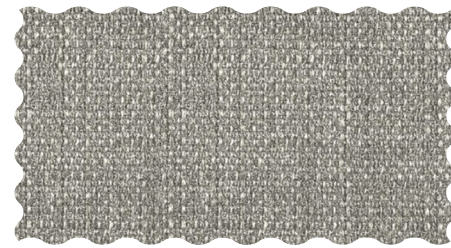

Silver

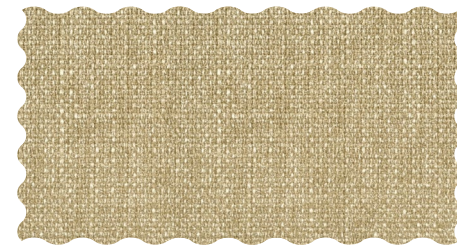

Hessian




Steel




Dove




Chablis

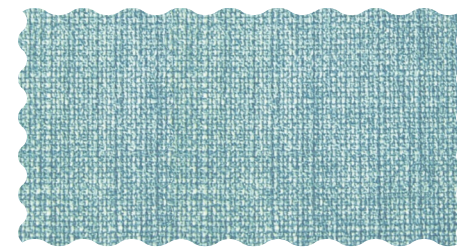


Sky

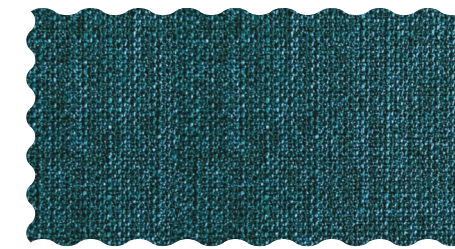


Ocean

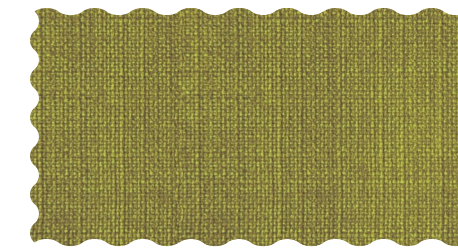


Fern green

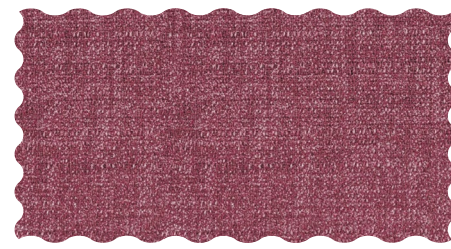


Lilac

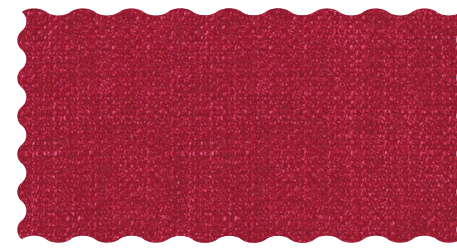


Cyclamen

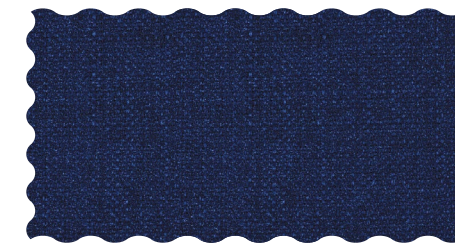


Royal

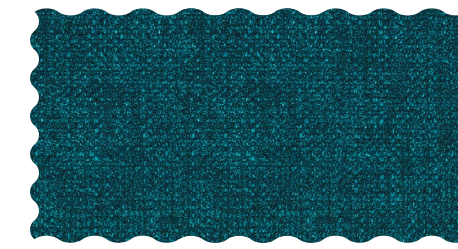


Teal

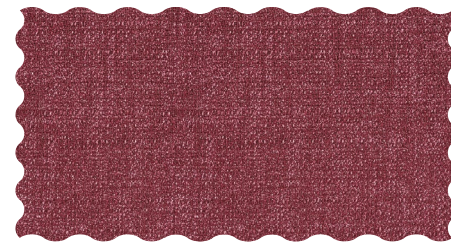

Peony

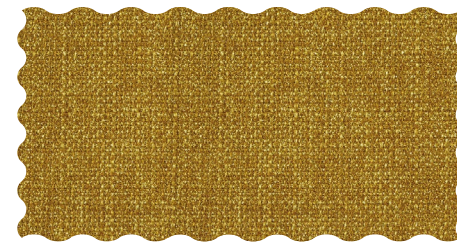

Gold

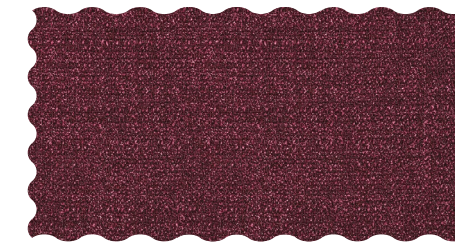

Fig



Mushroom

This is just a selection of our upholstery options, with many more available to choose from, including wipeable fabrics. Please ask our product specialist for details.

Accessories



Lifting pole



1050 mm width extension



Bed lever



High safety mat

Accessories



Folding side rails lowered




Folding side rails raised




Folding side rails with bumpers



Fabric side rails lowered

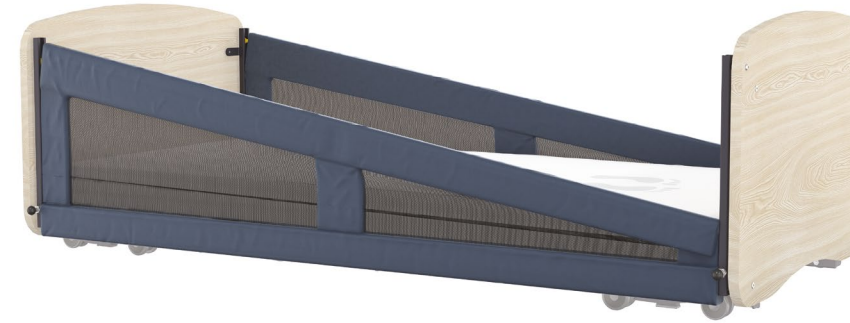

Fabric side rails partially raised

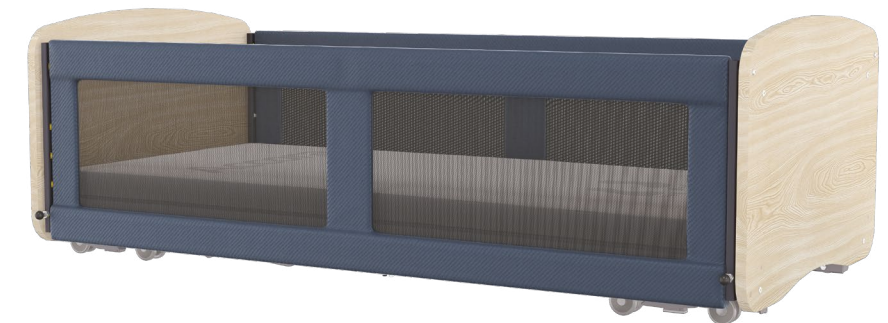

Fabric side rails raised

Technical specification

Overall dimensions (Arc headboard)	1015 mm width x 2090 mm length
Mattress size standard width Mattress size with width extension Mattress size with length extension	900 mm width x 2000 mm length 1050 mm width x 2000 mm length 2200 mm length
Mattress platform height	100 mm to 800 mm
UniSafe braking system	2 sets of brakes, each secure 4 castors
8 castors	4 lockable castors; 4 castors to aid manoeuvrability
Castor diameter	63 mm
Maximum lifting pole load	75 kg (12 st)
Auto regression	BodyMove technology
Backrest max. angle	70 degrees
Legrest max. angle	36 degrees
Trendelenburg function	12 degrees
Liquid ingress protection	IPX6
Safe working load	275 kg (43 st)
Max. patient weight	240 kg (38 st)
Warranty	Frame: 10 years Electric/parts: 3 years

The Empresa bed is part of our
Eminence™ Collection.

THE
EMINENCE
— COLLECTION —

Tailored by you. Made to order by Accora.

1. Lutheran Home. FloorBed: resident falls analysis. Interim report. February 2019

2. McCarthy M. Falls are leading cause of injury deaths among older people, US study finds. BMJ. 2016; 354. <https://doi.org/10.1136/bmj.i5190>/accessed March 2021

3. Data on file. February 2021

4. Smith J. The role of interior design in the care environment. <https://www.homesmiths.co.uk/interior-design-in-the-care-environment/> accessed November 2020

UK:

Charter House, Barrington Road
Orwell, Cambridge SG8 5QP

T: +44 (0)1223 206 100

E: info@accora.care

W: www.accora.care

IE:

38 Main Street, Swords, Co. Dublin,
Ireland, K67 E0A2

T: +353 (0)1 695 0614

E: info@accora.care

W: www.accora.care

* and ™ are trademarks of Accora Limited. As our policy is one of continuous improvement, we reserve the right to modify designs without prior notice. Colours may vary. All information is correct at time of printing and is subject to change.

Accora

BRO-FB6-002UK Rev03