
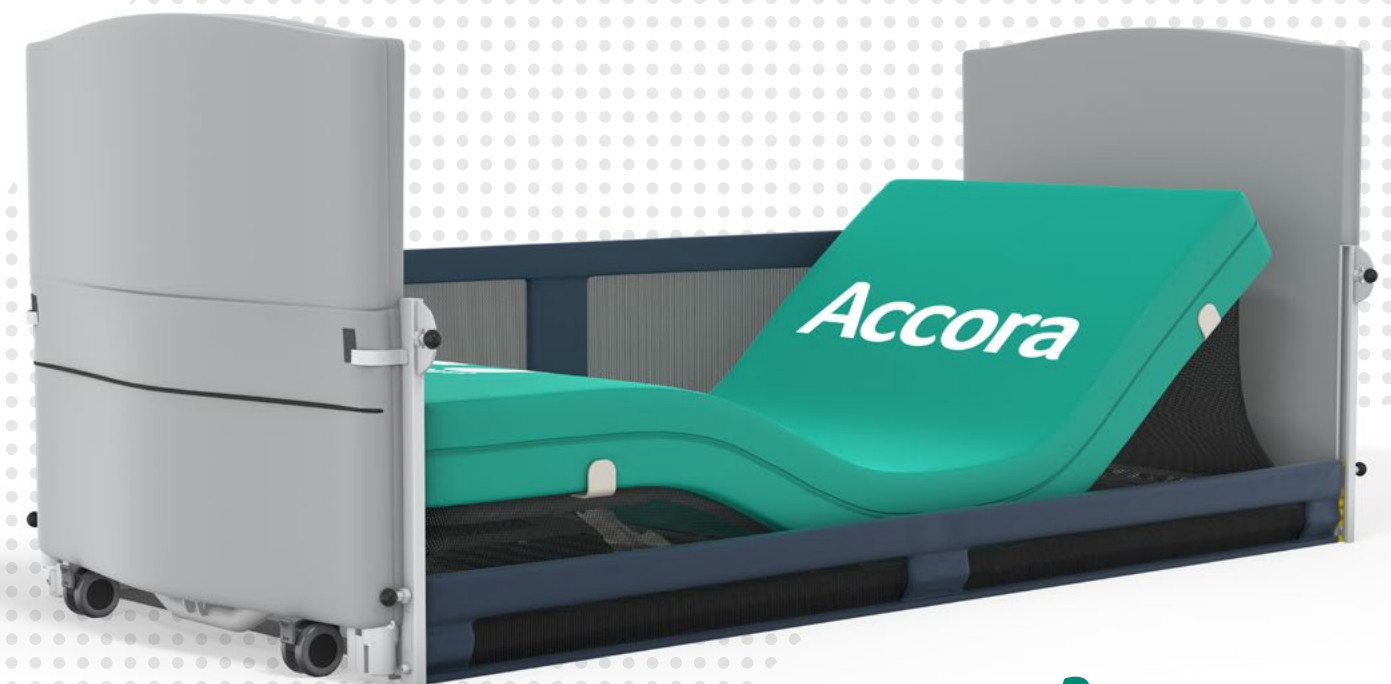


JuniorBeds

Flexible **solution** for a **safe**
sleeping **environment**

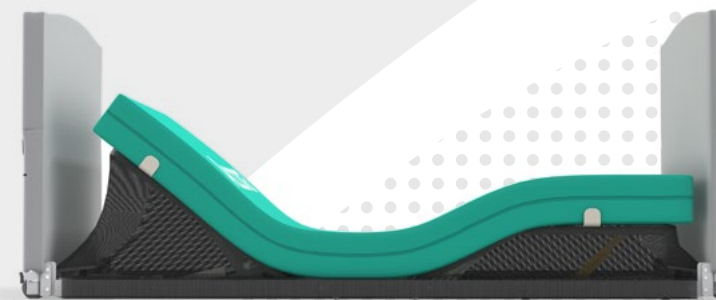

Accora

UK:

Accora Ltd

Charter House, Barrington Road, Orwell, Cambridge SG8 5QP
T: +44 (0)1223 206 100 **E:** info@accora.care **W:** www.accora.care

IE:

Accora Ltd

38 Main Street, Swords, Co. Dublin, Ireland, K67 E0A2
T: +353 (0)1 695 0614 **E:** info@accora.care **W:** accora.care

References

1. Data on file. 2. Zhu Hui-ping et al (2012) Disability, home physical environment and non-fatal Injuries among young children in China. PLOS ONE, 7(5): e37766. Available from <https://doi.org/10.1371/journal.pone.0037766> [Accessed 12th April 2021] 3. Hyndman B (26th March 2019) Preventing fall-related injuries in children: an environmental scan of resources and evidence-informed best practices. Ontario Neurotrauma Foundation. Available from <https://onf.org/ONFchildfallprevscanfinal2> [Accessed 12th April 2021]

* and ™ are trademarks of Accora Limited. As our policy is one of continuous improvement, we reserve the right to modify designs without prior notice.

BRO-FB4-001UK Rev02

Accora

JuniorBeds

The power of 3

Accora design is built around safety, practicality and affordability. **JuniorBeds** are a perfect blend of all three.



Suitable for use with **FloorBed 1**, **FloorBed 2** and **CommunityBed**, the **JuniorKit** suite of combined accessories provide a complete solution to meet the latest children's bed standard BS EN 50637:2017.

Creating a safer environment

These easy-to-fit accessories help transform an adult medical bed into a safe sleeping space for younger people and small adults.



FloorBed 1



FloorBed 2



CommunityBed


Flexible assets

Having the flexibility to make your standard Accora beds suitable for children and individuals with a smaller stature enables optimal use of your core stock and negates the need for special orders, whilst meeting your regulatory obligations.

It allows best use of your budget and time, too. And when the **JuniorKit** is no longer required, it can be recycled to benefit other service users.

Padded head and footboard bumpers help reduce impact and minimise the risk of injury.


Cost effective, simple and safe, the Accora bed range and **JuniorKit** are perfect for your service needs and your budget.

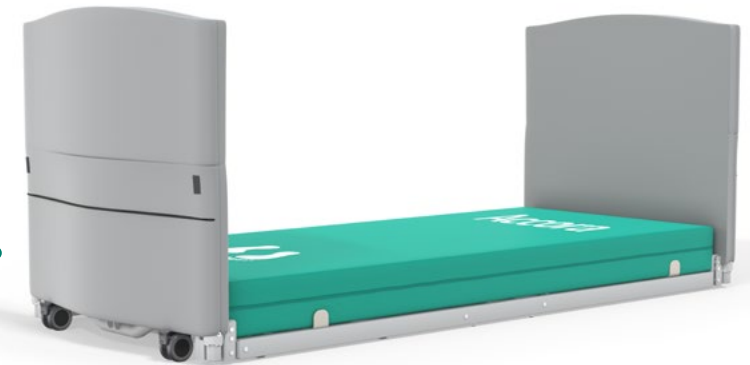

Built-in, lockable pocket for storing the handset to prevent tampering.



Safety sleeve encases the platform mechanism to help protect against risk of entrapment.

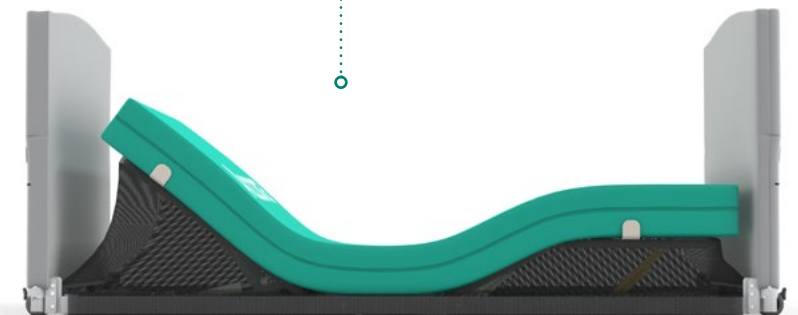

JuniorBeds

Design for living

Our industry-first **FloorBed** technology is all about changing lives.¹ Falls are common in children, but even more so in those with a disability.^{2,3}

Accora's **JuniorBeds** not only help make a young person's environment safer by reducing the likelihood of falls, they help them make the most of their abilities, encouraging independence and self-esteem.

And being able to adjust the height of the bed and its profile at the touch of a button helps take some of the strain out of caregiving, making carers' lives safer too.




Accora

JuniorBeds

Continuing support

The decision to use a specialist bed can be difficult at any time, but especially so with children. As a younger person outgrows their **JuniorBed**, the choice of accessories with the **FloorBed** and **CommunityBed** means they can be adapted to suit changing needs, without the disruption of having to provide a new bed.

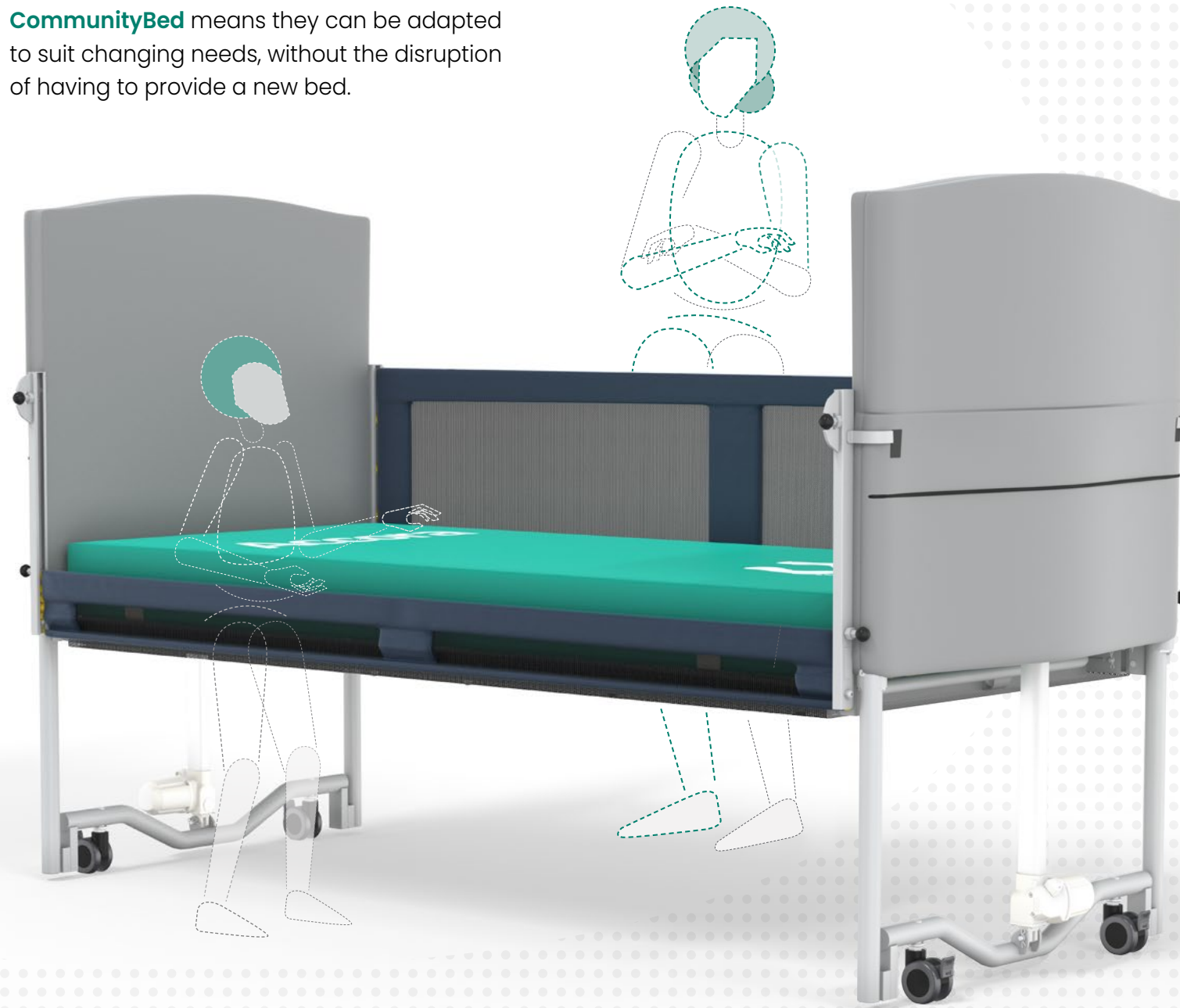

Choosing the right bed

Following these steps here will guide your choice of bed.

Is your service user someone who is more than 146 cm tall?
Weighs more than 40 kg?
Has a BMI of more than 17?

Yes

Consider standard Accora bed that meets BS EN 60601-2-52



FloorBed 1



FloorBed 2



CommunityBed

No

Is your service user less than 75 cm tall?

Yes

Consider a crib or alternative product

No

Does your service user require high side rails (over 45 cm high from mattress platform)?

Yes

Consider a cot or alternative product






No

Consider **FloorBed** or **CommunityBed** with **JuniorKit** to meet children's standard BS EN 50637:2017

JuniorBeds

How the FloorBed measures up

Comparing FloorBed 2
with a traditional nursing bed

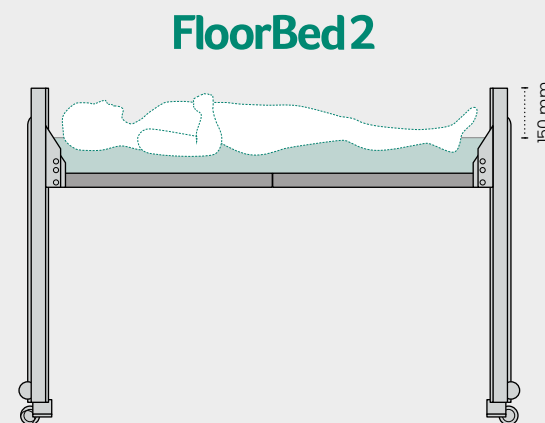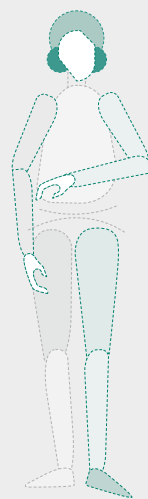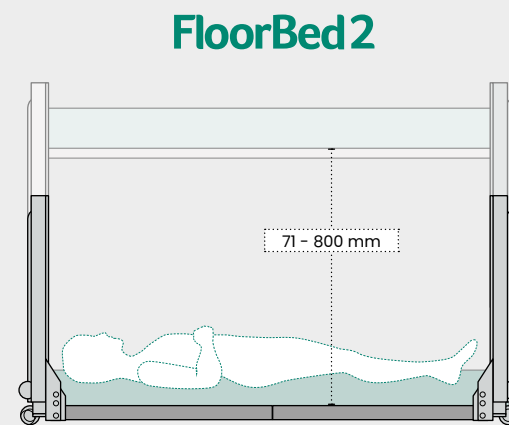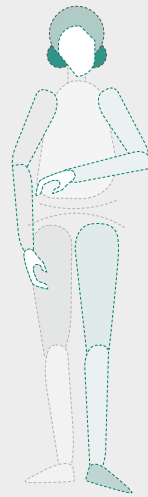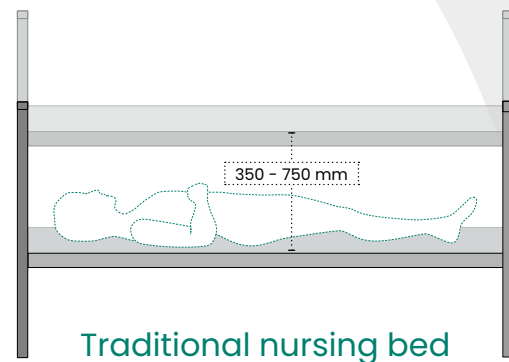

FloorBed and CommunityBed

Tailored for care – all the essential profiling functions and features to help assist the safe care and independence of children and small adults.

FloorBed 1

A true floor bed at just 71 mm high and offering superb user friendliness and durability, this market-leading bed has become a staple of equipment stores.


Its design makes this a good general-purpose bed, and for more vulnerable users, paediatric occupational therapists may prefer the benefit afforded by its higher head and footboards.



FloorBed 2

The same great floor level feature as **FloorBed 1**, but with the capability to be raised even higher to 800 mm where needed, to facilitate caregiving.

Also comes with low head and footboards for easier access for care tasks and to sliding sheets, allowing a single carer to reposition a person up or down the bed.

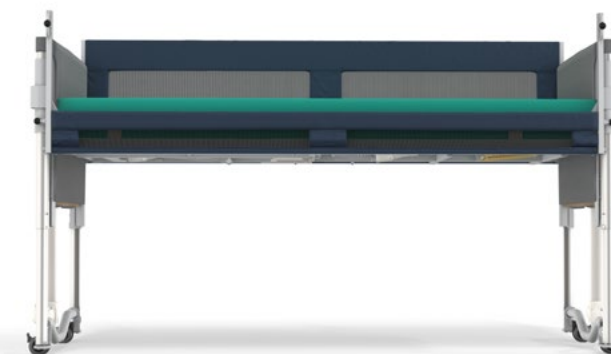

CommunityBed

A versatile “two-in-one” bed, which combines the 800 mm lifting height of a conventional profiling bed and the convenient height of a standard low bed, plus a maximum user weight of 180 kg.

Enables you to simplify your range and reduce stock holding, whilst still catering for the needs of smaller service users.




JuniorBeds

Bed accessories

Only these accessories are compatible for use with JuniorBeds.

Fabric side rails

With easy lock and lower mechanism. Fabric mesh allows visibility from within and outside the bed.

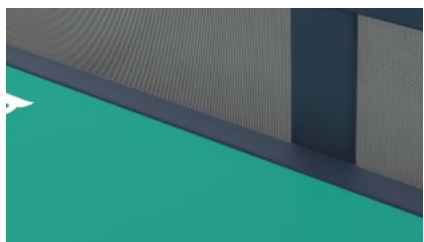

Integrated bolster

Forms in-built mattress infill. Side rails lower to below mattress height to facilitate transfers.

Firm edge mattress

Designed with a firm edge to help further reduce the risk of entrapment.



Easier transfers

Also helps service users orientate themselves in bed and can make transfers easier.

High safety mat

Helps reduce the risk of injury from bed falls for someone who is prone to falling.



Extended surface

When used with the bed at floor level, the high safety mat provides an extended surface to safely roll onto.

Technical specification



	FloorBed 1	FloorBed 2	CommunityBed
Overall dimensions	925 mm width x 2255 mm length (36.4 in x 88.8 in)	935 mm width x 2285 mm length (36.8 in x 90 in)	930 mm width x 2250 mm length (36 in x 89 in)
Mattress size	900 mm width x 2000 mm length (35.4 in x 78.7 in)	900 mm width x 2000 mm length (35.4 in x 78.7 in)	900 mm width x 2000 mm length (36 in x 79 in)
Safe working load	185 kg (29 st)	185 kg (29 st)	215 kg (34 st)
Service user weight	150 kg (24 st)	150 kg (24 st)	180 kg (28 st)
Mattress platform height	71 mm – 650 mm (2.8 in – 25.6 in) (excl. mattress)	71 mm – 800 mm (2.8 in – 31.5 in) (excl. mattress)	220 mm – 800 mm (8.6 in – 31.5 in) (excl. mattress)
Profiling functions	Backrest with auto-regression, kneebreak, Anti-Trendelenburg, Trendelenburg	Backrest with auto-regression, kneebreak, Anti-Trendelenburg, Trendelenburg	Backrest with auto-regression, kneebreak, Anti-Trendelenburg, Trendelenburg
Minimum height of service user	146 cm (4.8 ft) With JuniorKit 75 cm (2.5 ft)	146 cm (4.8 ft) With JuniorKit 75 cm (2.5 ft)	146 cm (4.8 ft) With JuniorKit 75 cm (2.5 ft)
Part number	NSB-0-FL1-200	NSB-0-FL2-200	NSB-0-CBI-100
Junior Kit			
Part number	JUNKIT-0-FL1-000	JUNKIT-0-FL2-000	JUNKIT-0-CBI-000
Comprises	Padded head and footboard bumpers Lockable pocket for handset Safety sleeve	Padded head and footboard bumpers Lockable pocket for handset Safety sleeve Protector plate	Padded head and footboard bumpers Lockable pocket for handset Safety sleeve
Note: JuniorBeds can only meet the new children's standard when installation of the full JuniorKit is complete, as specified here			